


NEW!
The best
just got better!

Geo-Mattress® Max

Premium Pressure Relieving Foam Mattress

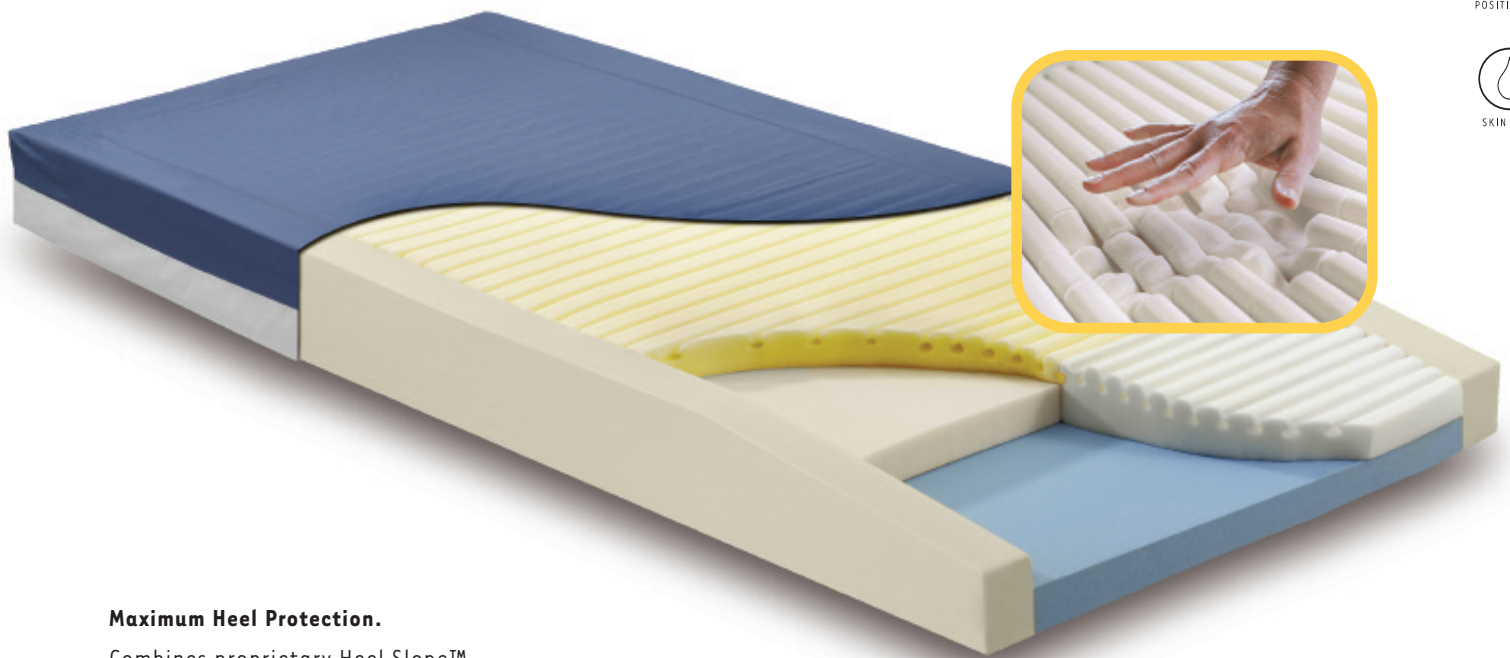

Maximum Heel Protection.

Combines proprietary Heel Slope™ and all-new viscoelastic “memory” foam to deliver unsurpassed offloading of vulnerable heels.

Maximum Safety.

All-new firm perimeter bolsters facilitate safer patient transfers, help eliminate dangerous gaps, and add stability for edge-of-bed sitting.

Maximum Performance.

Clinically proven Geo-Matt® anti-shear surface, multi-tiered construction, 7-year warranty and 400 lb. capacity make it our finest offering.



SURFACES



SEATING



POSITIONERS




SKIN CARE

SPANAmerica
Innovative Solutions.

Geo-Mattress® Max

Premium Therapeutic Foam Mattress

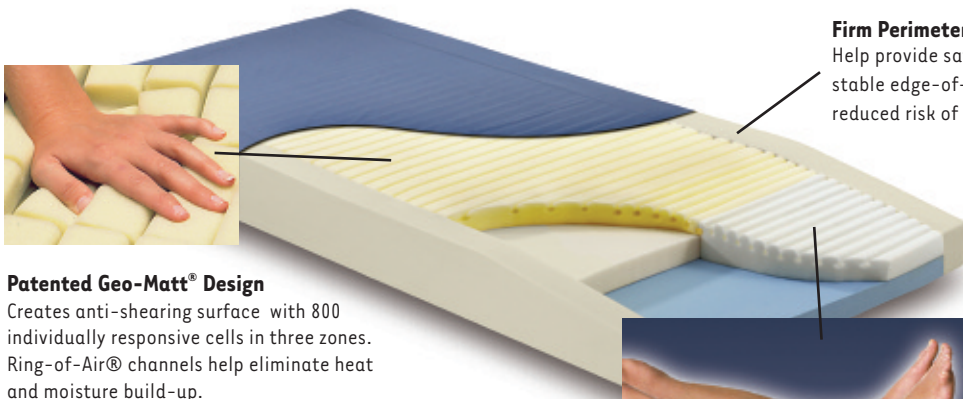

Firm Perimeter Bolsters

Help provide safer transfers, more stable edge-of-bed sitting, and reduced risk of gaps at siderails.

Patented Geo-Matt® Design

Creates anti-shearing surface with 800 individually responsive cells in three zones. Ring-of-Air® channels help eliminate heat and moisture build-up.

- **Multi-tiered, progressive resistance core.**
Premium, high-density foams are combustion-modified for conformance to Cal TB 117.
- **Fluid-repellent, flame-resistant top fabric.**
Anti-microbial and anti-fungal. Neoprene/butyl coated underside yields extremely low moisture vapor transmission rate (MVTR). Conforms to NFPA 702 Class 1 and Cal TB 117E.
- **Exclusive hidden zipper cover design.**
Flapless construction creates no collection points for contaminants. Combines with stain-resistant surface coating to ensure easy, thorough cleaning.
- **Flammability compliant.**
Standard versions conform to requirements of Cal TB 117 and 16 CFR Part 1632.4. Available in inner fire barrier versions for conformance to NFPA 101 Life Safety Code, ASTM E 1590, Cal TB 129 and 603, and Boston (BFD) IX-11. †
- **7-year, non-prorated warranty.**
Supports up to 400 lbs.
All materials latex-free.



By Design: Heel Protection Times Two.

Two distinct features to protect vulnerable heels, the site of roughly one third of all pressure ulcers¹:

Exclusive Heel Slope

Subtle taper redistributes load onto more pressure-tolerant lower legs while providing complete foot support. Documented to deliver heel pressures an average of 27% lower than leading competitors².

Viscoelastic foot section.


Luxurious, "memory" foam with individually-segmented Geo-Matt design. Enhanced immersion gently cradles heels.

1. Blaszczyk J. Make a Difference: Standardize Your Heel Care Practice. *Ostomy/Wound Management*. 1998;44(5):32-40.
2. McGrew, G. A Comparative Study of Interface Pressures on Two Hospital Polyurethane Foam Mattresses. *Clinical Reports: Series on Skin and Wound Care Management* 1997.


A Geo-Mattress for every need:




Geo-Mattress with Wings®
Patented, raised perimeter design.

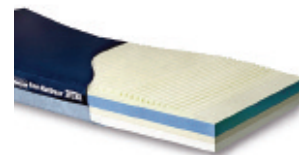



Geo-Mattress Atlas®
Bariatric support up to 600 lbs.

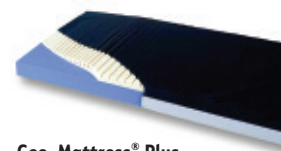

Geo-Mattress® Plus
Dual layer, 6-year warranty.

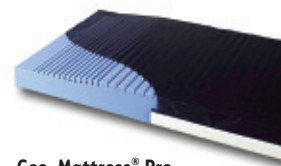

Geo-Mattress® Pro
Single layer, 5-year warranty.

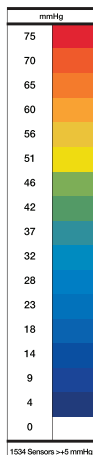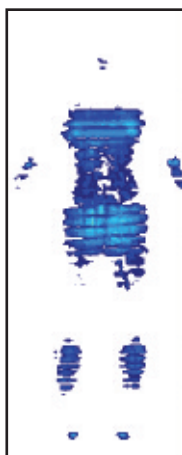

Computerized Pressure Mapping confirms exceptional pressure distribution, extremely low peak pressures.

Subject cited:
Male, 5'10", 175 lbs.
Supine, with pillow.
(Complete results available upon request).

Institutional Sizes		Consumer Sizes	
Item #	Description	Item #	Description
MX8435-13	84"L x 35"W x 6"H	* MX7539-13	75"L x 39"W x 6"H Twin
MX8035-13	80"L x 35"W x 6"H	* MX7554-13	75"L x 54"W x 6"H Full
* MX8032-13	80"L x 32"W x 6"H	* MX8060-13	80"L x 60"W x 6"H Queen
* MX7835-13	78"L x 35"W x 6"H	* MX8076-13	80"L x 76"W x 6"H King
MX7535-13	75"L x 35"W x 6"H	* allow additional lead time for these sizes	
Weight limit: 400 lbs.		Mattress weight: approx. 24 lbs.	
Warranty: 7 years.			

† For inner fire barrier versions, replace "-13" suffix with "-29".

SPANAmerica
Innovative Solutions.